

assisted bathing solutions

CI/SIB	(74.26)	X	(U35)

The Chiltern Invadex Shower Trolleys and Stretchers



Bringing care to the community

Shower stretcher and shower trolley



The Chiltern Invadex range of shower trolleys and stretchers are designed to allow easy washing of patients in the prone position.

Mobile and easy to manoeuvre, the shower trolley range makes transferring a user to the bathing area simple.

The Chiltern Invadex trolley combines performance and functional design to ensure the user is showered, dried and clothed without the need to lift manually, in line with manual handling regulations.

The shower/changing stretcher range (placed over a bath or in a shower area) features a convenient method of storage and is designed for ease of use by nursing staff.



T1 Mobile Shower Trolley

Shower trolley

The easy steer mechanism makes moving between rooms or hospital wards easy, with three locking castors to prevent movement when a patient is being transferred to the trolley. The fourth castor has directional locking to assist manoeuvrability.

Two convenient drop sides make transfer to and from the trolley a simple operation. The easy to use height adjustment mechanism allows a carer to move the trolley up and down to a convenient height, avoiding the need for stretching or leaning and can be raised or lowered from either side.

Safety

The safe and comfortable design features safety catches on the drop sides, lockable castors and a padded liner with head support.

Full size model

T1 Shower Trolley (1930mm x 750mm) will comfortably accommodate a large adult. The T2 is a fixed height trolley. The T4 is electrically operated by hand control.

Smaller top model

T3 Shower Trolley (1600mm x 700mm) with smaller top is available where space is limited.



T13 Shower bathing/changing unit
Tip-up wall mounted model



The Chiltern Invadex Shower/Changing Stretcher range is safe and sturdy providing reassurance for the bather and choice of height helps to make showering simple for the carer.

complete access to the bath or other facilities when required. A version with back support (T12) is available, which can be adjusted to seven positions, allowing the user to sit up if required.




Folds neatly away

For ease of storage the stretcher folds away neatly against the wall, allowing

Shower/changing stretcher



T12 Wall mounted tip-up shower/changing stretcher with adjustable back

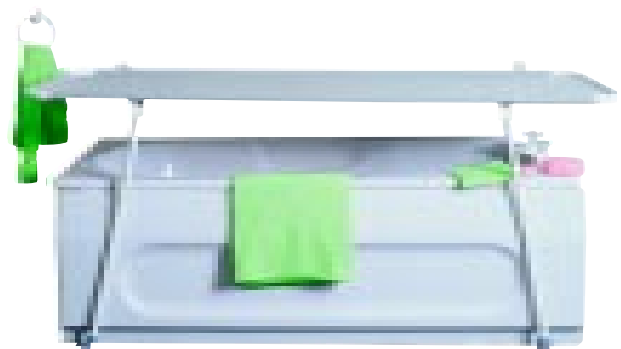

T11 Wall mounted tip-up shower/changing stretcher

Shower Trolley

- **Height adjustable**
Easy to use hydraulic height adjustment.
- **Convenience**
The trolley can be raised or lowered from either side, to a comfortable height for the carer with convenient drop sides for easier patient transfer.
- **Safety**
Safety catch fitted to drop sides. Directional locking on fourth castor to assist manoeuvrability.
- **Designed for the heavier patient**
Carries up to 190kgs (30 stone).
- **Easy steering**
Straight line steering device, allows easy manoeuvring from one area to another.
- **Braked castors**
Castors can be locked in position for added safety.
- **Comfort**
Padded bathing tray liner complete with head support for added comfort.
- **Easy drain**
Water drains away easily through an outlet hose.
- **No manual handling**
Users can be showered and dried without the need to lift manually.

Shower/Changing Stretcher

- **Model options**
The tip-up stretcher model is adjustable and can be locked in various positions to suit the user.
- **Easy store**
Folds neatly away against the wall when not in use.
- **Durable materials**
Cover made from tough pvc. Framework coated to prevent corrosion.
- **Easy drain**
Drain holes allow water to run away quickly and easily.
- **Made to order**
Special sizes can be made to order. Please ask for details.

accessories

Grab rails



Tip up shower seats




the shower trolley and stretcher range technical specification

Handing of T12 Shower Stretcher


The adjustable backrest on the T12 is on the carers right hand side as standard. Opposite hand made to order.

To arrange for Chiltern invadex to fit your selected stretcher please contact your nearest Regional Office.


Shower/Changing Stretcher Wall mounted tip-up models

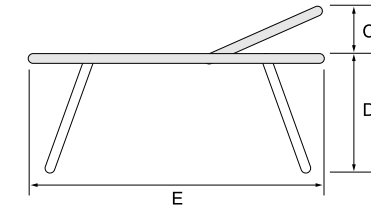

T12 Wall mounted tip-up with adjustable back

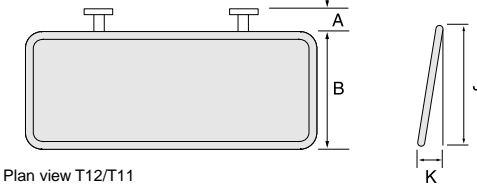


Plan view T12/T11

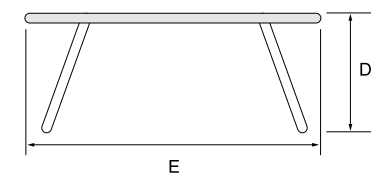

T11 Wall mounted tip-up

Shower Trolley


Mobile and Wall mounted tip-up models

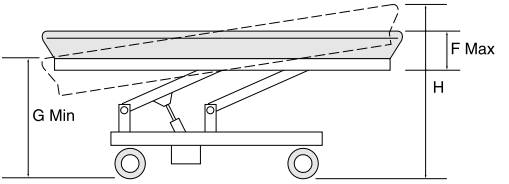

T1 Wheeled shower trolley



Plan view T1/T13

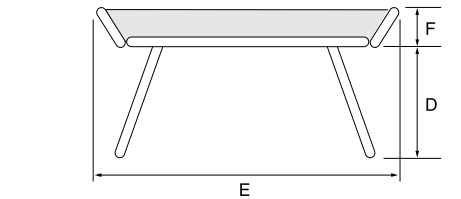

T13 Wall mounted tip-up shower bathing/changing unit

Technical Information - dimensions shown in mm and are for guidance only

	A	B	C	D	E	F	G	H	J	K
T1 Trolley	-	750	-	-	1930	200	580	960	-	-
T2 Trolley	-	750	-	-	1930	200	Fixed	960	-	-
T3 Trolley	-	700	-	-	1600	200	580	960	-	-
T4 Trolley	-	750	-	-	1930	200	580	960	-	-
T11 Stretcher	125	700	-	900	1900	-	-	-	700	110
T12 Stretcher	125	700	ADJ	900	1900	-	-	-	700	110
T13 Stretcher	175	750	-	900	1930	200	-	-	750	185

Quality

Chiltern Invadex is a service oriented company with more than 25 years experience in the caring market.

We design, manufacture, supply, install and service from sites throughout the United Kingdom and Ireland.

Chiltern Invadex products are also available through a network of overseas distributors.



FM 40009
FM 21085

Bicester and Dartford are BS EN ISO 9001 approved (FM 40009).

Bangor Northern Ireland is BS EN ISO 9002 approved (FM 21085).

Brochure Hotline

FREEPHONE 0800 917 0399


To request brochures on other products in our range

Chiltern Invadex Ltd
Chiltern House
6 Wedgwood Road Bicester
Oxon OX6 7UL
Tel: 01869 246470
Fax: 01869 247214

Overseas enquiries
Tel: +44 (0)1869 365545
Fax: +44 (0)1869 365567

Chiltern Invadex (NI) Ltd
5 Balloo Way
Balloo Industrial Estate Bangor
County Down BT19 7QZ
Tel: 02891 450497
Fax: 02891 270046

Chiltern Invadex Ltd
66-68 Manchester Road
Chapel-en-le-Frith High Peak
SK23 9TP
Tel: 01298 816366
Fax: 01298 816106

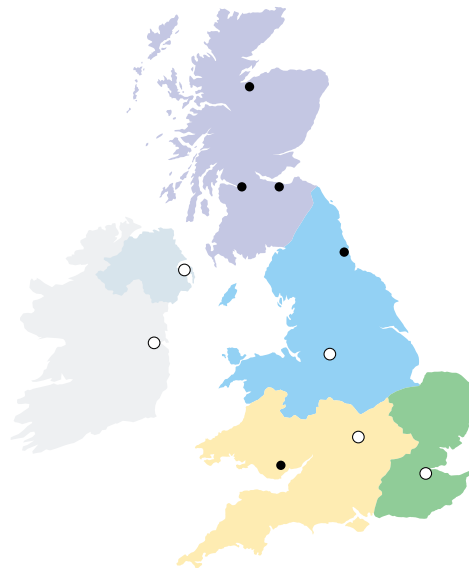

Chiltern Invadex Ltd
Units 5&6 Orbital One
Green Street Green Road
Dartford Kent DA1 1QG
Tel: 01322 284623
Fax: 01322 284579

Chiltern Invadex Ltd
St. Kierans Enterprise Centre
Furze Road Sandyford Industrial Estate
Foxrock Dublin 18
Tel: 01 286 2726
Mob: 086 829 7541

Web: www.chilterninvadex.co.uk

As part of our commitment to new product development we reserve the right to alter specifications without prior notice.

Published by : Chiltern Invadex Ltd, 6 Wedgwood Road, Bicester, Oxon, OX6 7UL Registered in England No: 1182024 ISSUE 1 AUG 1999
In Northern Ireland : Chiltern Invadex (NI) Ltd, 5 Balloo Way, Balloo Industrial Estate, Bangor County Down, BT19 7QZ Registered No: NI 19958

