




SHOWER TOILET CUBICLE RANGE

SHOWER CUBICLE SELECTOR - FRONT ENTRY TOILET CUBICLES


Use this shower tray selector to choose the cubicle base size required.

FRONT ENTRY 100MK2 SHOWER TOILET CUBICLE

Standard doors

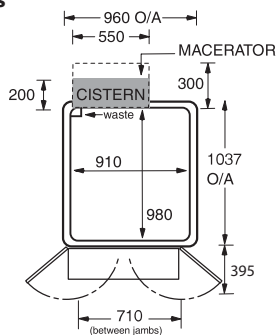

Bifold doors

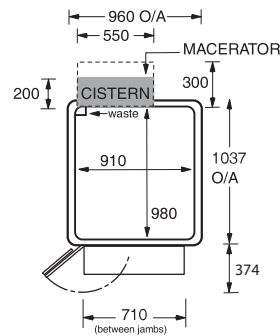


Ramp size
950mm Width x 280mm Depth


Note:
Left hand bifold illustrated. Right hand version available.

FRONT ENTRY SHOWER TOILET CUBICLE - 101

Standard doors

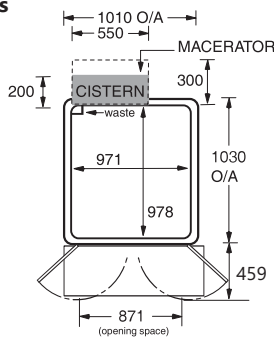

Bifold doors

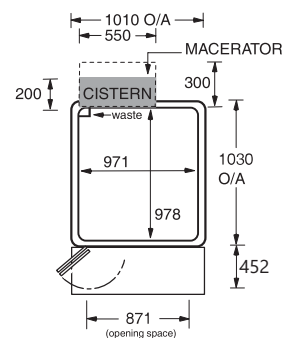


Ramp size
1010mm Width x 430mm Depth


Note:
Left hand bifold illustrated. Right hand version available.

FRONT ENTRY SHOWER TOILET CUBICLE - 120

Standard doors

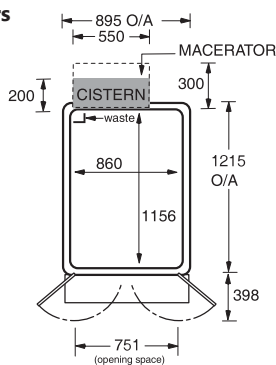

Bifold doors

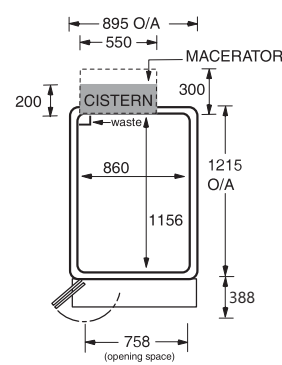


Ramp size
890mm Width x 240mm Depth


Note:
Left hand bifold illustrated. Right hand version available.

FRONT ENTRY SHOWER TOILET CUBICLE - 120PLUS

Standard doors

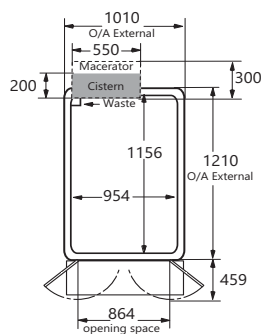

Bifold doors

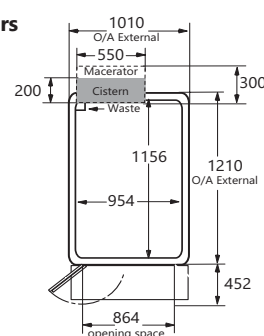

Ramp size
1010mm Width x 470mm Depth


Note:
Left hand bifold illustrated. Right hand version available.

SHOWER CUBICLE SELECTOR - SIDE ENTRY TOILET CUBICLES


Use this shower tray selector to choose the cubicle base size required.

SIDE ENTRY MK2 SHOWER TOILET CUBICLE

Standard doors

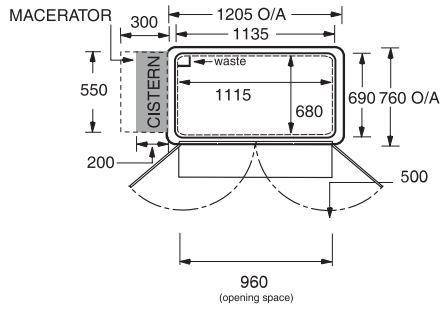

Bifold doors

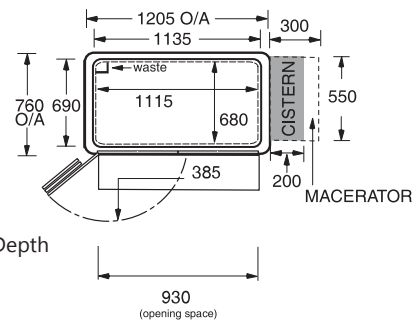

Ramp size
1160mm Width x 250mm Depth

Note:
Left hand bifold illustrated. Right hand version available.

- Notes:**
1. Add 50mm at rear of cubicle to allow for waste hose - pumped waste
 2. Add 100mm to width to allow for shower unit services
 3. Allow as much space as possible to gain access to services for routine maintenance.



SHOWER CHAIR CUBICLE TOILET COMPATIBILITY



Cubicle Model	Shower Chair Model											
	A02	A06	A08	A11	A13	A20	A21	TA18AP	TA18SP	TA21AP	TA21SP	TiS
100	✓	✓	✓	✗	✗	✗	✗	✓	✗	✓	✗	✗
101	✓	✓	✓	✗	✗	✗	✗	✓	✗	✓	✗	✗
120	✓	✓	✓*	✓	✗	✗	✗	✓	✓	✓	✗	✗
120Plus	✓	✓	✓	✓	✗	✗	✗	✓	✓	✓	✗	✓
Side Entry	✓	✗	✗	✗	✗	✗	✗	✗	✗	✗	✗	✗

* 21" Model - standard door option only
TiS model cannot be tilted when in the cubicle





HOW TO ORDER

Cubicle kit contents:



- White shower tray with ramp
- Full height enclosure pack (half height carer assist panels available on request)
- WC pan and cistern
- Plastic basin (fillable via shower hose)
- Door set - half height or bifold doors
- Curtain track, hooks & white shower curtain
- 2x Grab rails
- Fitting instructions

Extras

Ceramic basin & lever taps shown on p25
Shower unit shown on p25
Macerator shown on p26

	Option 1	Option 2	Option 3
Cubicle description	Half height door	Bifold (Left hand)	Bifold (Right hand)
100MK2 Cubicle Toilet with pumped waste	Part no: SHTCU100PW		
100MK2 Cubicle Toilet with gravity waste	Part no: SHTCU100G		


Technical information

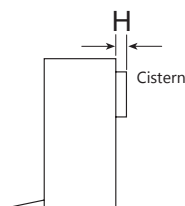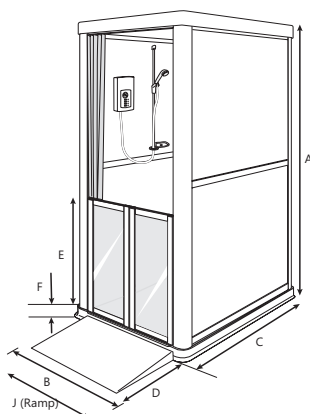



Ref	Dimensions
A	2030
B	950*+
C	1037
D	280
E	810
F	50
H	200
J	950

*Add 100mm if side controls fitted.
+Add 75mm if basin fitted.

Standard doors

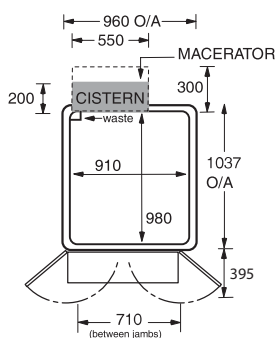

Bifold doors

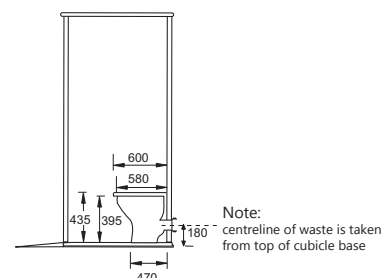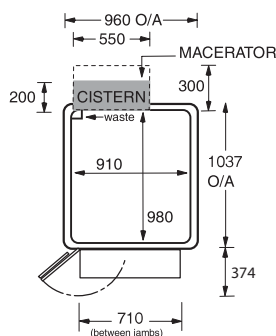

Notes: Dimensions are for doors fully open.

Internal dimensions: The width is measured from acrylic panel to acrylic panel. The depth is taken from the back acrylic panel to the inner face of the door.



Shower chair compatibility

See page 19 for shower chair compatibility options





HOW TO ORDER

Cubicle kit contents:



- White shower tray with ramp
- Full height enclosure pack (half height carer assist panels available on request)
- WC pan and cistern
- Plastic basin (fillable via shower hose)
- Door set - half height or bifold doors
- Curtain track, hooks & white shower curtain
- 2x Grab rails
- Fitting instructions

Extras

- Ceramic basin & lever taps shown on p25
- Shower unit shown on p25
- Macerator shown on p26




	Option 1	Option 2	Option 3
Cubicle description	Half height door	Bifold (Left hand)	Bifold (Right hand)
101 Cubicle Toilet with pumped waste	Part no: SHTCU101PW		
101 Cubicle Toilet with gravity waste	Part no: SHTCU101G		

Technical information

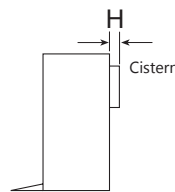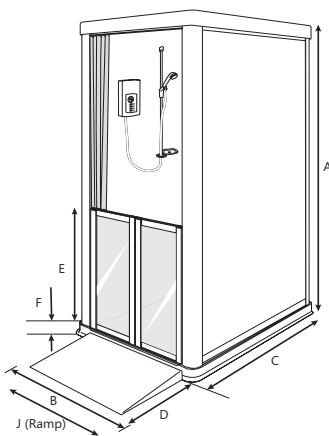

Ref	Dimensions
A	1900
B	1010*+
C	1030
D	430
E	810
F	50
H	200
J	1010

*Add 100mm if side controls fitted.
+Add 75mm if basin fitted.

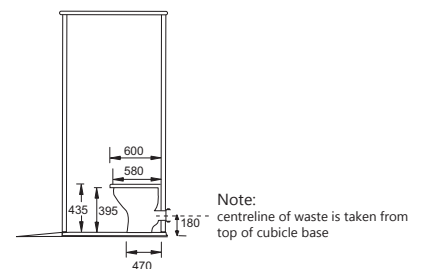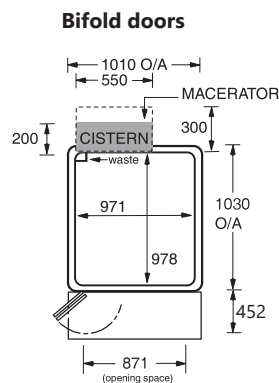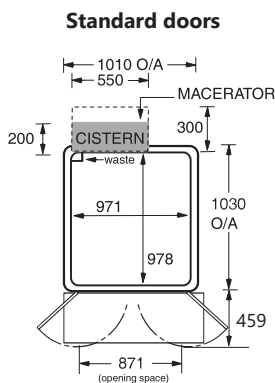

Notes: Dimensions are for doors fully open.

Internal dimensions: The width is measured from acrylic panel to acrylic panel. The depth is taken from the back acrylic panel to the inner face of the door.



Shower chair compatibility

See page 19 for shower chair compatibility options





HOW TO ORDER

Cubicle kit contents:



- White shower tray with ramp
- Full height enclosure pack (half height carer assist panels available on request)
- WC pan and cistern
- Plastic basin (fillable via shower hose)
- Door set - half height or bifold doors
- Curtain track, hooks & white shower curtain
- 2x Grab rails
- Fitting instructions

Extras

Ceramic basin & lever taps shown on p25
Shower unit shown on p25
Macerator shown on p26

	Option 1	Option 2	Option 3
Cubicle description	Half height door	Bifold (Left hand)	Bifold (Right hand)
120 Cubicle Toilet with pumped waste	Part no: SHTCU120PW		
120 Cubicle Toilet with gravity waste	Part no: SHTCU120G		

Technical information

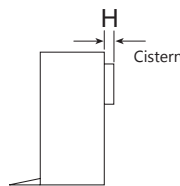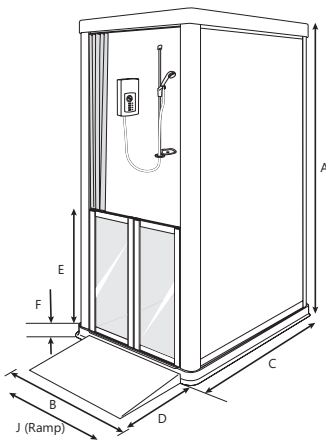


Ref	Dimensions
A	1900
B	890*+
C	1215
D	240
E	810
F	50
H	200
J	890



Notes:

*Add 100mm if side controls fitted.
+Add 75mm if basin fitted.

Standard doors

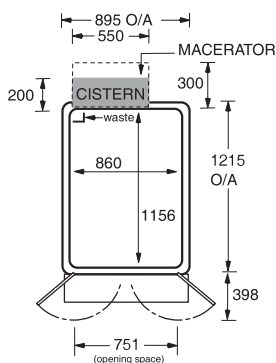

Bifold doors

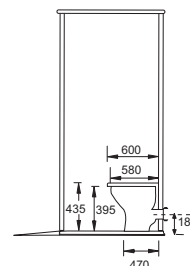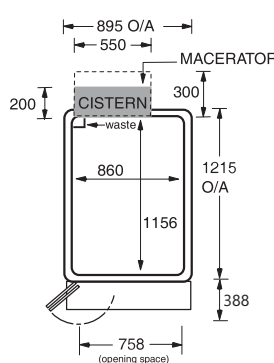

Note: centreline of waste is taken from top of cubicle base

Notes: Dimensions are for doors fully open.

Internal dimensions: The width is measured from acrylic panel to acrylic panel. The depth is taken from the back acrylic panel to the inner face of the door.



Shower chair compatibility

See page 19 for shower chair compatibility options





HOW TO ORDER

Cubicle kit contents:



- White shower tray with ramp
- Full height enclosure pack (half height carer assist panels available on request)
- WC pan and cistern
- Plastic basin (fillable via shower hose)
- Door set - half height or bifold doors
- Curtain track, hooks & white shower curtain
- 2x Grab rails
- Fitting instructions

Extras

Ceramic basin & lever taps shown on p25
Shower unit shown on p25
Macerator shown on p26

	Option 1	Option 2	Option 3
Cubicle description	Half height door	Bifold (Left hand)	Bifold (Right hand)
120Plus Cubicle Toilet with pumped waste	Part no: SHTCU120PWPLUS		
120Plus Cubicle Toilet with gravity waste	Part no: SHTCU120GPLUS		

Technical information

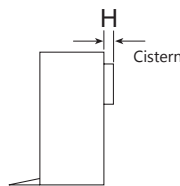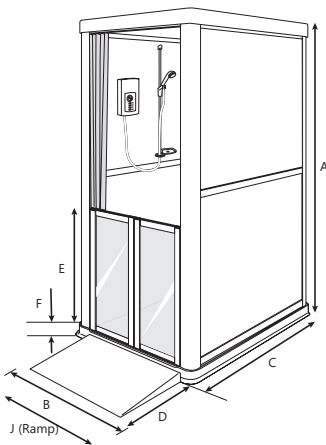


Ref	Dimensions
A	1900
B	1010
C	1215
D	470
E	810
F	50
H	200
J	1010



Notes:

*Add 100mm if side controls fitted.
+Add 75mm if basin fitted.

Standard doors

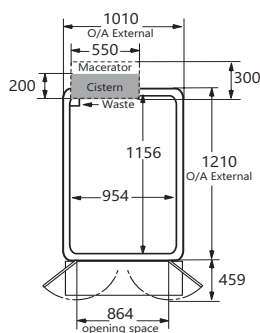

Bifold doors

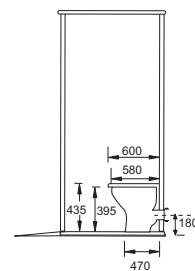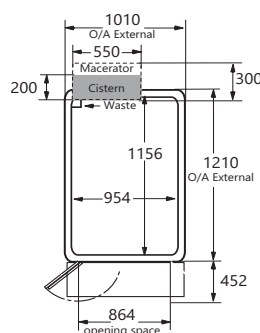

Note: centreline of waste is taken from top of cubicle base

Notes: Dimensions are for doors fully open.

Internal dimensions: The width is measured from acrylic panel to acrylic panel. The depth is taken from the back acrylic panel to the inner face of the door.



Shower chair compatibility

See page 19 for shower chair compatibility options





HOW TO ORDER

Cubicle kit contents:



- White shower tray with ramp
- Full height enclosure pack (half height carer assist panels available on request)
- WC pan and cistern
- Plastic basin (fillable via shower hose)
- Door set - half height or bifold doors
- Curtain track, hooks & white shower curtain
- 2x Grab rails
- Fitting instructions

Extras

Ceramic basin & lever taps shown on p25
Shower unit shown on p25
Macerator shown on p26

	Option 1	Option 2	Option 3
Cubicle description	Half height door	Bifold (Left hand)	Bifold (Right hand)
Side Entry MK2 Cubicle with pumped waste	Part no: SHTCUSE2PW		
Side Entry MK2 Cubicle with gravity waste	Part no: SHTCUSE2G		

Technical information

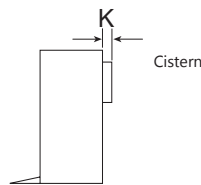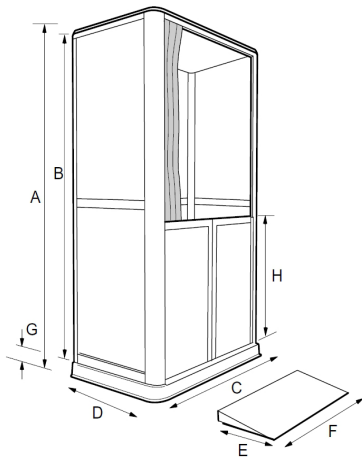


Ref	Dimensions
A	2030
B	1950
C	1160*+
D	740
E	250
F	1160
H	810
K	200



Notes:

*Add 100mm if side controls fitted.
+Add 75mm if basin fitted.

Standard doors

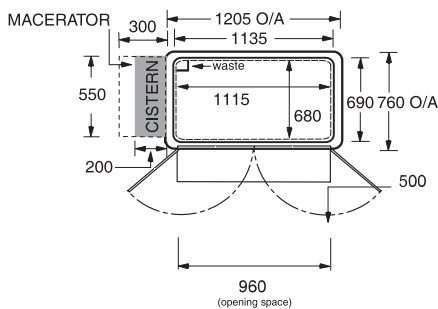

Bifold doors

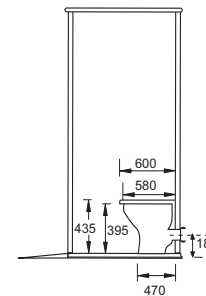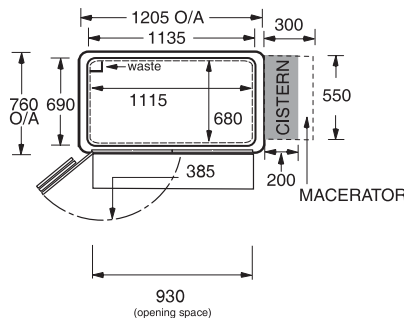

Note: centreline of waste is taken from top of cubicle base

Notes: Dimensions are for doors fully open.

Internal dimensions: The width is measured from acrylic panel to acrylic panel. The depth is taken from the back acrylic panel to the inner face of the door.



Shower chair compatibility
See page 19 for shower chair compatibility



CUBICLE EXTRAS



Ceramic basin and lever taps are available as an alternative to the standard plastic basin.

This option must be specified at time of order to enable reinforcement of the side panel.

Cannot be retro-fitted

code

SHSU029

product

Ceramic cubicle basin & lever taps



Mira Advance ATL Flex Extra Electric Shower

- Mira Advance ATL Flex 9kW electric shower offers a flexible solution for those who are looking for safety and increased mobility within the showering area
- Extended temperature control lever, 2m hose and additional handset holder for flexible showering heights.
- Compatible with the Whale Instant Match Pump see page 39
- 2 Year guarantee

code

IDWMA05A

product

Mira Advance ATL Flex - Electric 9 Kw shower heater, incl. 1 metre riser rail and 2 metre hose, separate wall mounted handset holder, standard riser rail, and easy touch temperature control lever



Mira Excel Thermostatic Mixer shower

- Showering with precise thermostatic control, recognising and compensating for changes in water temperature to provide an extremely safe and satisfying shower
- Easy installation features save time and installation costs
- Includes maximum temperature stop for safe showering
- 5 year guarantee

code

IDWME01A

product

Mira Excel - Thermostatic Mixing Valve (handset, riser rail and hose included)

SANITARY SYSTEMS



Saniplus small bore sanitary system

- Breaks down toilet waste from a shower, wash basin or W.C for pumping away through a small bore (22mm) pipe
- A powerful motor enables the unit to pump the waste and water up to 5m vertically, 100m horizontally or a lesser combination of both, where conventional large bore pipework is impractical
- Features a new case with sound reduction and integrated carbon filter
- Can be used with Whale pump option

code	product
SH0411	Saniplus sanitary system

Technical information

- Motor induction motor (no carbon brushes or gears to wear out). Tested over 30,000 cycles.
- Watts (Consumption) 400w
- Electrical Supply 220-240V 50Hz
- Flow rate at max head l/min 40 Horizontal Discharge 100 metres
- Vertical Discharge 5 metres
- Gravity fall on horizontal installations 1:100
- Discharge Pipe Diameter 22mm
- Materials Polypropylene, C.R.P., stainless steel, Neoprene Type of WC Pan (not included)
- Short spigot 100mm HO outlet (BS EN 33/37) Motor RPM 2800 RPM
- Normal Running Time WC 10-20 secs approx or as required, to discharge bath, shower etc
- Dimensions (vertical height is measured from floor level) 512W x 305 O/H x 200D (requires 210mm behind spigot)

Designed to fit BS EN 33/37 WC pans.

As close coupled cistern dimensions vary, it may be necessary, with certain designs, to utilise a spacer between cistern and wall


Saniplus small bore sanitary system (shown with Whale pump option)

Saniplus small bore sanitary system with Whale pump


Where the situation demands, or if the pipe run is greater than that specified for the Whale pump a Saniplus small bore sanitary system may be used. The discharge from the Whale pump must go into the top of the Saniplus small bore sanitary system.

Do not use the bottom entries.

It is preferable to have two separate discharge lines to waste and drainage when used in the same room as a Saniplus small bore sanitary system. Any failure of the Saniplus small bore sanitary system will not be detected by the Whale pump.

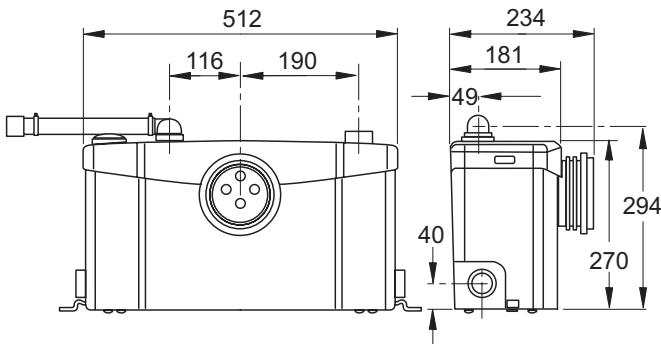

Pumping distance illustration

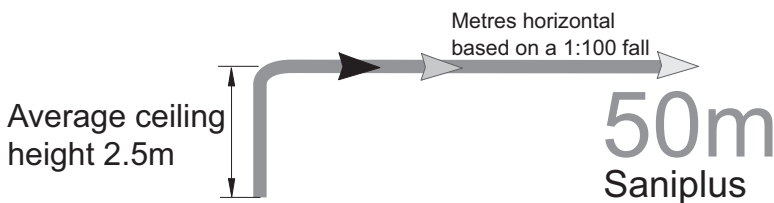

NB: Vertical lift must precede horizontal run

Note: If you increase 'horizontal' pipe fall to full 1:40 gravity fall, then the horizontal distance can be increased to virtually any distance required (also increase pipe size)

Macerator shown in relation to a Chiltern Invadex Shower toilet cubicle

